BébéSounds®

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bebesounds.com

©2004 Listro Associates

*Please Read The Following
To Enjoy Our Product!*

- ***Please have patience. This product requires a little time but the results are wonderful.***
- ***What you will hear through the BébéSounds Monitor will be the NATURAL sounds of your baby.*** Your doctor uses Doppler ultrasound equipment that produces artificial noise. Our BébéSounds Fetal Monitor amplifies the actual and natural sounds of your baby.
- Use a good battery. An old one will interfere with your BébéSounds Fetal Monitor. We recommend a new 9-volt alkaline battery.
- The best time to listen is 3-4 hours after eating and when you are relaxed.
- The Monitor must be used in a very quiet room.
- The Monitor must be placed on bare skin. Turn the Monitor on, place the listening cone against your bare stomach (start two inches below and several inches to the right of your belly button), and press the activation button to listen. Move the Monitor an inch at a time but be sure not to press the activation button as you move it.
- Placing the volume control at a low to moderate level will produce the best results and avoid static.

- You might not hear sounds immediately so don't be disappointed. You must have patience when trying to hear your baby's sounds. The more you use the BébéSounds Monitor, the better you will become at finding the sounds. The baby's position will also greatly affect your success. The best position is with his/her back to your belly button. Your success will also depend on how relaxed you are. You may have to try several times before you will hear your baby's sounds, particularly the heartbeat, which may be heard as early as the fifth month. Other sounds are heard earlier.
- Your baby's heartbeat does not sound like your own. A fetal heart beats from 120-180 times per minute. It will sound like a soft, fast beating drum, or a galloping horse.
- You will also hear kicks, hiccups, and nutrients passing through the placenta. Kicks sound like uneven thumps, hiccups like two rapid drumbeats, and the nutrients like a whooshing sound.

BébéSounds is a Personal Care Product
and cannot be returned to the retailer.

*All Exchange And Refund Requests Should
Be Directed To Our Help Desk*

**If you are having difficulty, have questions, or
need assistance, call our Help Desk toll free at
1-888-232-6476 Mon-Fri from 9:00AM-5:00PM EST.**

Please Read On

INTRODUCTION


Please read the following material carefully to maximize your enjoyment of your BébéSounds Prenatal Heart Listener.

It is important to understand that many factors affect what you will hear with your Fetal Monitor. The position of your baby, the baby's weight, your weight, the stage of pregnancy, and where you place the BébéSounds Monitor, all help to determine what you will hear.




For example, when listening for the fetal heartbeat, the sound will be loudest in the third trimester and if your baby is positioned with his/her back to your belly button. Your baby may naturally shift position, thereby creating the situation where you will hear the baby one moment but not the next. However, once you've passed your second trimester, your baby will be growing very rapidly and you will hear the heartbeat and all other sounds on a much more regular basis.

Don't worry! It's quite normal that there will be times you will not be able to hear your baby's sounds due to your baby's position or stage of development. This is especially true of the heartbeat which is highly dependent on your baby's position and is a very soft, fast, beating sound like a drum or galloping horse.

KIT COMPONENTS

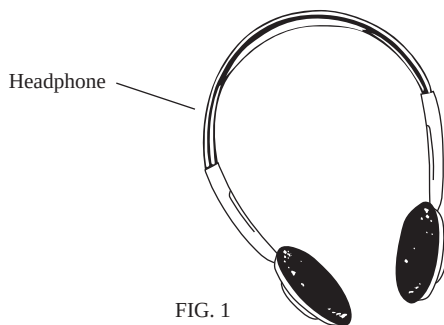Instruction Manual

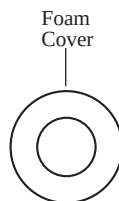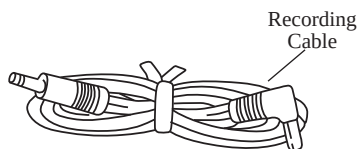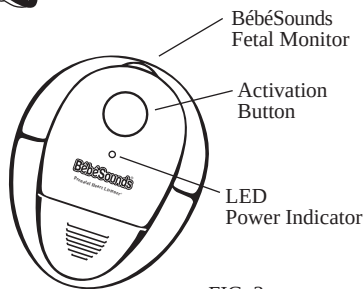IDENTIFYING THE FETAL SOUNDS YOU WILL HEAR

You will hear kicks, hiccups, nutrients passing through the placenta, and, of course, the heartbeat. Each one is distinctive and being able to identify them will significantly enhance your experience.

We think the best way to describe each sound is as follows:

Kicks:	A series of uneven thumps
Hiccups:	Two rapid, soft drumbeats
Nutrients:	A whooshing sound
Heartbeat:	A very soft rapidly beating drum or a galloping horse

A fetal heart beats from 120-180 times a minute while the heart of a pregnant woman beats approximately 100 times a minute. Knowing this should help prevent you from confusing your own heartbeat with your baby's.

WHEN YOU WILL HEAR THESE SOUNDS

The heartbeat may be heard as early as your fifth month. You will need to listen very carefully for your baby's heartbeat as it may be extremely difficult to hear it over the beating of your own heart. In the earlier months of your pregnancy, the sound of your own bowels may also interfere with your ability to hear the heartbeat. While BébéSounds is specially designed to tune out these sounds, they can still get in the way.

That is why the best time to listen is 3-4 hours after a meal and when you are relaxed and have a slower heartbeat.

You may be able to hear the passing of nutrients through the placenta, kicks, and hiccups as early as your second trimester. Remember, as your baby grows and gets stronger, so do the sounds he/she creates. You are also becoming more familiar with the sounds so you will be able to identify them more easily.

THE NATURAL SOUNDS OF THE BÉBÉSOUNDS FETAL MONITOR VS. THE ELECTRONIC NOISE CREATED BY A DOCTOR'S DOPPLER

You will notice that the sounds you hear through the BébéSounds Monitor are different from those you have heard at your doctor's office. This is because your doctor probably uses Doppler ultrasound equipment that creates electronic noise based on sound waves. BébéSounds, however, is designed to generate the sounds naturally by using the absolutely safe Bell stethoscope method of amplification. Again, the heartbeat heard through the Monitor sounds more like the soft, fast, beating of a drum, while the Doppler sound is more of a gushing noise. Please keep in mind that it will take a little practice to pick up this faint sound so listen very carefully.

BébéSounds is a Personal Care Product
and cannot be returned to the retailer.

*All Exchange And Refund Requests Should
Be Directed To Our Help Desk*

**If you are having difficulty, have questions, or
need assistance, call our Help Desk toll free at
1-888-232-6476 Mon-Fri from 9:00AM-5:00PM EST.**

RECORD YOUR BABY'S SOUNDS AND YOUR OWN HEARTBEAT

With our recording cable and your own tape recorder you can record your baby's sounds. We also encourage you to record your own heartbeat prior to your baby being born. The sound of your heartbeat will have a calming and soothing effect on your newborn baby.

OPERATION OF BÉBÉSOUNDS PRENATAL HEART LISTENER (Refer to Figures 6 and 7)

1. Make sure the foam cover is in place on the listening cone.
2. Remove battery compartment door and install a new 9 volt battery (not included) keeping the ribbon under the battery and exposed for easy removal. Press battery forward and down. We recommend a new alkaline battery as an old one will severely interfere with the performance of your Monitor.
3. Insert the headphone plug into either jack on the Monitor.
4. Turn on BébéSounds using the ON/OFF volume control switch. The red LED will light to indicate the unit is in the ON position.
5. Use the BébéSounds Monitor on bare skin in a quiet room.
6. Hold the conical-shaped detector to your stomach or back and press and hold down the activation button. Adjust the volume to optimize the listening level and amplification. Start with a low to moderate setting to avoid static.
7. Release the activation button and gently move the Monitor one inch at a time to different places until you find the best position for listening to your baby's sounds. Start two inches below and several inches to the right of your belly button.

TOP VIEW AND BOTTOM VIEW GRAPHIC

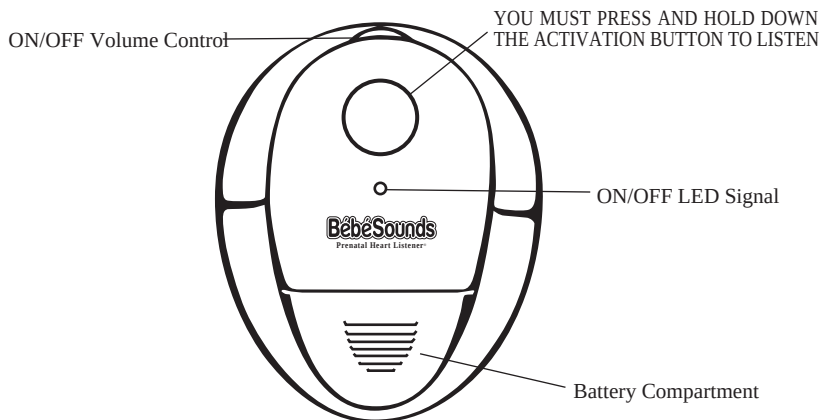
FIG. 6

FRONT VIEW

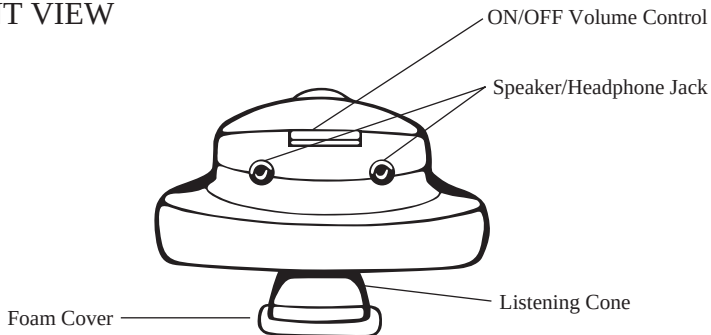
FIG. 7

A FEW TIPS ABOUT USING BÉBÉSOUNDS FETAL MONITOR

- The BébéSounds Monitor has powerful amplification capabilities required to pick up the quiet sounds produced by your baby. To protect your ears and to prevent feedback (squealing sounds), DO NOT press the activation button unless the listening cone is pressed against your body.
- Increasing the distance between the headphones and the Monitor may help you avoid interference or feedback (squealing sounds).
- There are many different sounds being produced by your baby and your body. Increasing volume also increases these sounds. Keeping the volume at a moderate level usually produces the best results, therefore, adjust the volume carefully.
- Your baby's heartbeat is a very low volume sound and is very difficult to detect until late in the second trimester. Patience and very careful listening are required to locate the best position to hear the heartbeat. It is perfectly normal to hear it on some days and not others due to your baby's position and movements. It is also normal that you will hear it one moment and not the next. Don't worry! It is very natural for your baby to shift its position.

TO RECORD FROM THE BÉBÉSOUNDS MONITOR

1. Plug the L-shaped end of the recording cable (see Fig. 4) into either jack of the BébéSounds Fetal Monitor and the straight end into the line input of your tape recorder. (Some portable recorders do not provide this input.)
2. Position the Monitor on your bare skin so that you can hear the sounds you wish to record.
3. Adjust the recording level on your tape recorder if necessary.
4. Hold the Monitor very still and begin recording.

Remember, the sound of your heartbeat has a very calming and soothing effect on your newborn baby. You may want to record your own heartbeat before your baby is born. You will then be able to play back your recorded heartbeat to your baby after birth. We recommend a 30-minute recording. (You will have to reduce the volume when you record your own heartbeat since it is much louder than your unborn baby's.)

The Prenatal Heart Listener[®] is a personal care product and cannot be returned to the retailer.

*All Exchange And Refund Requests Should
Be Directed To Our Help Desk*

**If you are having difficulty, have questions, or
need assistance, call our Help Desk toll free at
1-888-232-6476 Mon-Fri from 9:00AM-5:00PM EST.**

RECORD AND EMAIL YOUR BABY'S HEARTBEAT WITH YOUR COMPUTER

1. Locate the recording cable. Plug the L shaped end of the recording cable into one of the available jacks on the Prenatal Heart Listener® monitor. Plug the other end into the microphone jack on your computer. This is indicated by a picture of a microphone near the jack or the jack is color coded pink on the back of your computer.
2. On your computer go to START, PROGRAMS, ACCESSORIES, ENTERTAINMENT (or MULTIMEDIA), SOUND RECORDER.
3. Set up your Prenatal Heart Listener® for listening to your baby. Press the red recording button in the SOUND RECORDER program to begin recording a 2 minute sound file. Press the Prenatal Heart Listener® monitor against your tummy and hold down the activation button on the monitor. You are now recording your baby's sounds.
4. After 2 minutes the recording will finish. Name and save your sound file to an appropriate folder. You can now email your sound file to family and friends.

TROUBLESHOOTING

If unit is not working check the following:

1. Activation button not held down during listening.
2. Volume control turned down.
3. Battery and/or headphone connection not complete.
4. Foreign matter in listening cone.
5. Weak battery (ON/OFF LED signal will dim).
6. Feedback/squealing noises – turn down volume and/or move headphones away from BébéSounds. Press activation button only when BébéSounds Monitor is against your body.
7. Foam cover not in place on listening cone.

DISCLAIMER BÉBÉSOUNDS PRENATAL HEART LISTENER IS NOT A MEDICAL DEVICE

BébéSounds is not a medical device, and should not be used or relied upon in any way for diagnostic or other medical purposes. Any concerns about your pregnancy should be addressed to a medical practitioner. **BébéSounds Prenatal Heart Listener is not intended to be and should not be used as a substitute for regular prenatal care by an appropriate medical practitioner.** You are encouraged to seek such prenatal care in connection with your pregnancy. The BébéSounds Prenatal Heart Listener should be used only in accordance with the instructions provided herein, and the manufacturers and distributors of BébéSounds Prenatal Heart Listener are not responsible for inappropriate use. The manufacturer expressly disclaims liability for incidental, consequential, or special damages of any nature.

LIMITED PRODUCT WARRANTY

BébéSounds Prenatal Heart Listener, distributed by UNISAR Inc., has a limited warranty against defects for a period of 90 days from the original date of purchase.

IN-WARRANTY PROCEDURE

In the event service is required due to defect or malfunction during the warranty period, UNISAR Inc. will replace or repair the product under warranty. The repair or replacement of the defective item will be free of charge. Repair or replacement of the defective item is subject to verification of the malfunction or defect when delivered to:

Saro Import Export, Ltda.
General Diaz Polier, 89
28006 MADRID

Please enclose dated proof of purchase with your return.

The Prenatal Heart Listener[®] is a personal care product and cannot be returned to the retailer.

*All Exchange And Refund Requests Should
Be Directed To Our Help Desk*

**If you are having difficulty, have questions, or
need assistance, call our Help Desk toll free at
1-888-232-6476 Mon-Fri from 9:00AM-5:00PM EST.**

NOTE: The warrantor expressly disclaims liability for incidental, consequential, or special damages of any nature. (Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you.) This warranty does not cover (1) A defect that has been caused by repair or modification, (2) A defect caused by negligence, accident, or misuse of the product, (3) A defect caused by incorrect installation and/or defective batteries, or (4) Damage resulting from improper operation.

Any implied warranties shall be limited in duration to the terms of this warranty. (Some states do not allow limitations on how long an implied warranty lasts, so the above limitation may not apply to you.) This warranty gives you legal rights which may vary depending on which state you reside in.

NOTE: This warranty is valid only for products sold in the United States. If you have purchased your product outside of the United States, please see your local retailer or distributor for service.

IMPORTANT: Prior to returning BébéSounds Prenatal Heart Listener, replace the existing battery with a new one. Defective or exhausted batteries are the leading cause of problems with this item and are not covered by this limited warranty.

BébéSounds Prenatal Heart Listener is subject to United States patents and related foreign counterparts.

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Other Products by BébéSounds®

Angelcare® Movement Sensor w/ Sound Monitor
Stop worrying while your baby sleeps!™ Detects baby's slightest movements even while she sleeps. If absolutely no movement is detected for 20 seconds, an alarm sounds. The monitor alerts you that your baby has stopped moving and is absolutely still. Can also be used as just a sound monitor.
Suggested Retail: \$74.99

Color Portable Video & Sound Monitor
Our 2.4 GHz color LCD monitor has 4 channels and a 450 foot range. Both the camera and monitor are rechargeable. Features infrared night vision, voice activation, power saving mode, belt clip, and AV adapter for connecting your TV. *Suggested retail \$199.99*

Nursery Air Purifier
Helps remove odors from the nursery. Emits negative ions, which attach themselves to allergens, dust, airborne bacteria, and remove them from the air. Effective in a room up to 200 square feet. HRPF washable filter.
Suggested Retail: \$39.99

Call 1-800-430-0222

Other Products by BébéSounds®

Ear Thermometer

Our unique Ear Thermometer doesn't need a probe cover, takes just one second to get an accurate reading, and stores the last 10 readings in memory. Temperature measured in either Fahrenheit or Celcius. Includes storage case and battery. Suggested Retail: \$24.99

Portable MiniFridge

Keeps your food and beverages cold (41°F) or hot (149°F). Will hold tall or wide-necked baby bottles, baby food jars, sandwiches, snacks, or even six 12-ounce cans. Includes car adapter and home AC adapter.

Suggested Retail: \$59.99

Extra Headset for Use with FirstSounds® and BébéSounds® Prenatal Heart Listeners®
Enables a second person to listen to your unborn baby.

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**If you are having difficulty, have questions, or
need assistance, call our Help Desk toll free at
1-888-232-6476 Mon-Fri from 9:00AM-5:00PM EST.
Please ask for Millie Cupeles**

UNISAR Inc.
15 West 36th Street
New York, NY 10018
Tel: 1-888-232-6476